INN  
912

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Dental water jet by PASTA DEL CAPITANO produces short water steams so that it efficiently removes dental plaques and food leftovers from interdental spaces, reducing bacterial charge. Regular irrigation helps to prevent gum discomforts and improves dental and oral hygiene. Always use water jet after brushing teeth.

Read carefully manual instruction before using the appliance. Keep this manual for all life duration of appliance for future references

## **GENERAL INFORMATION AND SAFETY INSTRUCTIONS**

When using electric devices, especially with children, carefully read general safety instructions.

### **READ CAREFULLY THE FOLLOWING INSTRUCTIONS BEFORE USING DEVICE:**

- Before connecting device to plug, verify if voltage specified on the product label corresponds to the voltage of your home electric socket.
- Never submerge device into water or other liquids.
- Don not handle device with watered hands.
- Do not use the device next to bath tube or shower.
- Do not place the device in spots where can easily fall into the sink or shower.
- Only use provided accessories.

- Fill in water tank with water or fitting brawls for this kind of item.
- Do not use brawls with iodine tincture, bleach or oils.
- Do not put any object in the crevices.
- Do not use the device if it doesn't work correctly, damaged or fallen into water. In case of incorrect function address to an authorized client assistance centre.
- Before charging the device verify if USB cable is integrate and is not damaged
- Do not open, unhinge, or attempt to repair the device. Dismantle only in case of disposal of the appliance.
- Do not use the device outdoor, near a nebuliser or flammable substances.
- The device is designed for domestic use only.

## **FIRE HAZARD**

- The use of device can interfere with peacemaker or other medical equipment. Consult primary care physician before using water jet.
- In order to avoid bacteria growth into the tank or in the spout,

completely empty the tank after use to eliminate water from tube.

- Use the device following the instructions in this manual and directions of your dentist.
- Do not push the spout against gums or teeth.
- Do not direct water flow under the tongue, behind ears, into the nose or other sensible body areas. The water jet can cause redness and damages if headed to particularly sensible body areas. Carefully follow instructions, every different use from description does not conform to described use in the instruction manual and has to be considered inappropriate. The manufacturer declines every responsibility for damages by an improper use of the device.
- Keep the product out of the reach of children.
- This item can be used by children or people with reduces physical, mental and sensorial capacities, or without experience and knowledge only under supervision or if they received instruction about the safely use of devices or if they understand dangers associate to the use of device.
- Cleaning operations and maintenance should not be affected by unsupervised children.
- This device is not a toy. This product is not designed for being


- used by children without a direct supervision of an adult
- Packaging material (plastic bags, cartons, polystyrene ecc.) should not be left at the reach of children since they are potential sources of danger and should be disposed according to current laws.
  - Make sure that children do not play with the appliance.
  - This appliance is not intended for professional use in dental offices, use only to clean teeth and gums, in case of pain or bleeding contact the dentist. It is not a medical device.

### **TRADITIONAL ORAL HYGIENE**

Traditional toothbrush cannot clean away food particles and harmful bacteria.

The toothbrush could only touch a limited area under the gum line.
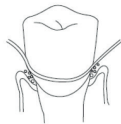
### **HYGIENE WITH INNOLIVING WATER JET**


Innoliving Dental water jet can clean away bacteria and food particles hiding under gums and inter-dental lines.
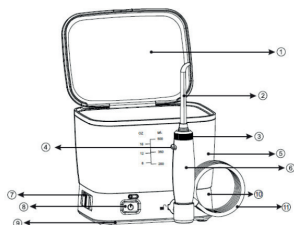
## PRODUCT DESCRIPTION

1. Water tank lid
2. Jet tip
3. Tip rotating knob
4. Jet tip eject button
5. Water tank
6. Handle
7. Water pressure switch
8. ON/OFF button
9. Base
10. Body
11. Water pipe

PART	DESCRIPTION	QUANTITY	IMAGE
1	Main unit	1	
2	Jet tips	2	
3	USB cable	1	
4	Handle	1	
5	Bag	1	

## INSTRUCTIONS

**WARNING:** Please fully fill in the Water Tank for the first use and completely let out from dispenser into the sink. Please pay attention to the following instructions:

### Fill in the tank


Remove the water tank and water tank lid from the base

Fill pure water into the water tank place it on the base and lightly press it down in order to assure it is correctly fixed.


### Inserting and removing the tip

The two jet tips are placed in the accessory compartment. Open the lid to take the jet tip and the handle.

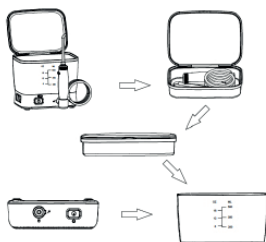
Insert the end of the handle until is blocked. To remove the jet lip push the lateral release button on the handle.

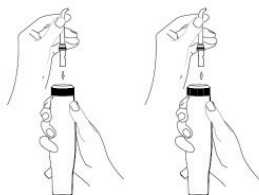
### Inserting handle on the base

Insert the handle on the base as shown in figure, then rotate handle 90° in clockwise to block it in upright position. Make sure the handle is strongly fixed. To remove the handle from the base, rotate the handle 90° in clockwise and disconnect from the base.

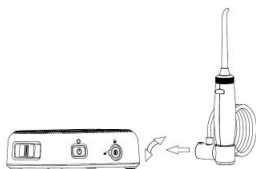
### Regulating water pressure

Move the ON/OFF button on the base side to increase or decrease water pressure.
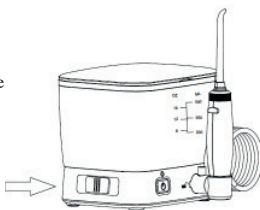
### How to use water jet


Lightly bow on the sink and insert the handle in oral cave, directing the jet tip towards the teeth. Push ON button on the base and ON/OFF button on the handle. Rotate ring on the handle to change the direction of the water jet.
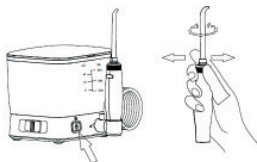
## Working time and battery

**Timer:** the oral irrigator has a preset working time of 150 seconds. For longer use, press the power button and ON/OFF switch on the handle.

The white light indicator starts to flash when there is low battery. It is recommended to charge the battery, only when it's low.  
To charge the product insert cable on an USB port, as shown in figure below.

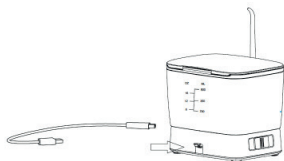
**When item is charging, white light will flash slowly. When light is steady battery is fully charged.**

**CAUTION: do not use the item, when it's charging the battery.**

## USING TIPS

For the first use of water jet it is suggested to use lower pressure level. Bow on the sink and direct water jet towards gums and teeth. Keep mouth slightly open to let water flow on the sink and to avoid water sprays to go in every direction. Direct water jet towards gums 'side, forming a 90° edge. Move long water jet on gums' side staying for some seconds on each tooth. For a better dental hygiene start from molars to incisors.

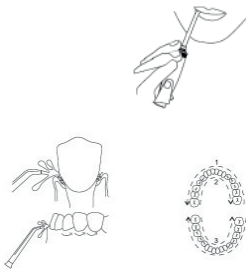Push ON/OFF button to switch off the

device before removing jet tip from mouth.

### After using

Empty laying water from tank. Switch on device for some seconds to clean away residual water on the tank and jet tip. Repeat this after every use to avoid bacteria growth. If it is no effected when device is not used for a period of time, bacteria can growth on water leftovers.

If your dentist suggests the use of a chlorhexidine solution, let water slide towards jet tip after use to avoid the substance to lay on the device.

### Cleaning

Switch off the device and remove water tank and jet tips. Clean the tank and jet tips with warm water. Clean the slot with a damp cloth. After cleaning rub with a dry cloth.

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# IPX7

Water resistant

This product complies with all the applicable European regulations

Protection class III

**INFORMATION TO USERS** according to Legislative Decree No. 49 of March 14, 2014 “2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE)“The crossed-out dustbin symbol indicates that the product at the end of its life must be

collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promote recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014.

Built-in batteries must be removed from products by recovery and recycling center. The batteries used in this device must be disposed of in special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps

preventing potentially negative consequences on the environment and human health.

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