INN  
010

**IDROPULSORE  
PORTATILE CORDLESS**

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**INSTRUCTION MANUAL**

*Thanks for purchasing INNOLIVING dental water jet Pasta del Capitano INN-010. This oral irrigator uses a penetrating water jet to clean the spaces between your teeth and to reach hard to clean areas that conventional toothbrushes may struggle to reach. It will also massage the gums.*

**IMPORTANT SAFETY INSTRUCTION**

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

**BEFORE USING THE APPLIANCE, READ THE FOLLOWING INSTRUCTIONS CAREFULLY:**

- Do not dismantle this product at any time other than when disposing of it.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Do not immerse in water or other liquid.
- Do not handle with wet hands.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Only use tips and accessories that are included.
- Fill reservoir with water.
- Do not drop or insert any foreign object into any opening on the product.
- Do not use this product if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. In case of doubts, please address to a specialised centre. Do not use mouthwashes/rinses that contain iodine, bleach, or oils.
- Do not use outdoors or operate where aerosol products are being used or oxygen is being used or administered. This is a FIRE HAZARD.
- Use of this product may interfere with the operation of pacemakers or other internal medical devices. If you use one of these devices you should consult your physician before using this irrigator.
- To avoid bacterial growth and contamination in water, always empty the water tank completely after use and clear the pump system by running the unit for a few seconds when the tank is empty.
- Never push the nozzle against teeth or gums
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage and keep the user manual for future reference.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this product or any other oral hygiene aid.
- Keep the product and its accessories out of reach of children. This product should only be used by children when supervised by adults.
- This product is not a toy. Children should be supervised to ensure that they do not play with the appliance.
- Never try to repair the product by yourself, otherwise the warranty will be not valid anymore.
- This product is not designed as a device for professional use in dental practices
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

**SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE**

## PRODUCT DESCRIPTION

- 1) Tip
- 2) Power switch/pressure control
- 3) Water
- 4) Accessory compartment
- 5) Battery compartment


## THE PACKAGE INCLUDES:

- 1 CORDLESS PORTABLE DENTAL WATER JET
- 2 INTERCHANGEABLE HEADS

## BEFORE USE

There may be moisture inside the main water unit on receipt of your Irrigator. This is moisture left from the distilled water used to quality control the item before it is sold. Insert 2xAA alkaline batteries into the battery compartment at the base of the unit. This compartment can be opened with a coin. Close the battery cover securely after installing the batteries.
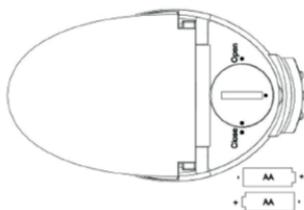
## GETTING STARTED

- 1) Insert the spray nozzle firmly to the top of the unit. The spray nozzle is stored in the front of the device. Lift up the lid, take the nozzle and insert it into the device.
- 2) Remove the blue tank and fill it with water. The water tank can be filled with cold or warm water (no hotter than 40°C). Reinsert the water tank into the device.
- 3) Insert the irrigation nozzle into your mouth, while you are leaning over a sink or basin. Hold the unit in an upright position with the nozzle pointing toward your teeth.
- 4) Gently push the switch upwards until it locks into the “High” or “Low” position. The pump will start operating. We recommend that you start off with the “Low” setting and move to the “High” setting as your gums get used to the massage motion and strengthen.
- 5) The nozzle can be rotated to reach all areas of your mouth. Your mouth should remain slightly opened while cleaning, allowing water to flow out.
- 6) Ensure that you keep the device in an upright position. If you turn the irrigator upside down, it will not take water from the tank.

**ATTENTION: fill the tank (3) with water, removing it from the main unit. Reposition the tank, when it's full.**

**ATTENTION: do not pour the water into the accessory compartment (4), to fill the water tank. Water leakage into inner components and battery room might occur in case of the water is dropped in that compartment, causing damages and malfunctioning.**

### Attention:

Never fill the tank with hot water. Gums can be very sensitive and easily damaged by hot water. Never turn the unit on before putting the nozzle in your mouth. Never push the tip of the nozzle against teeth or gums. Never use this device for cleaning any other part of the body such as ear or eyes. Please note that those with sensitive gums may experience light bleeding the first few times they use this irrigator. If bleeding persists, please consult your doctor.

**If water doesn't come out from the oral irrigator in case of first use, please suck the air from the tip to activate the water jet.**

### Using Mouthwash and Other Solutions

Oral Irrigator can be used to deliver mouthwash and antibacterial solutions. Ask your dental professional to recommend the appropriate solution for you.

Note: use of some antibacterial solutions could shorten the life of your oral irrigator. After using any special solution in the tank, take a moment after every use to rinse the unit. Partially fill the tank with warm tap water. Aim the tip into your sink. Turn on the unit and let it run until empty. This will help keep the tip from clogging and help your unit to last longer.

### AFTER USE

With the irrigation nozzle still attached, open the water tank by the top cap and drain away any remaining water. Switch on the unit for a few seconds to pump out any water that remains in the pump tank or nozzle. We recommend that you do this after every use to avoid encouraging bacteria. Remove the batteries from the unit, in case of long no-use.

### CLEANING AND MAINTENANCE

- 1) Switch off the unit, remove the blue water tank by sliding it downwards and remove the irrigation nozzle.
- 2) Wash the tank in warm water. Wipe the white housing with a damp cloth. Never use caustic or abrasive cleaning agents as they may damage the device.
- 3) Particular care should be taken when washing the hose connections in the tank.
- 4) Never submerge the unit.
- 5) After washing wipe with a dry cloth and store it in a cool and dry place, away from heat sources.

## TECHNICAL FEATURES

Code	Pasta del Capitano INN-010
Description	Dental water jet
Power	2xAA 1.5V alkaline batteries (included)

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Batterie: 2x 1.5V AA alkaline batt.

**MADE IN CHINA  
LOT.**

This product complies with all the applicable EU directives.

**INFORMATION TO USERS** according to Legislative Decree No. 49 of March 14, 2014 “2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE)”.

The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014.

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

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Prodotto sotto licenza ufficiale di Farmaceutici Dott. Ciccarelli S.p.A.