
innoliving

easy tech, easy life

**BILANCIA DIGITALE
PESAPERSONA
ULTRASLIM**

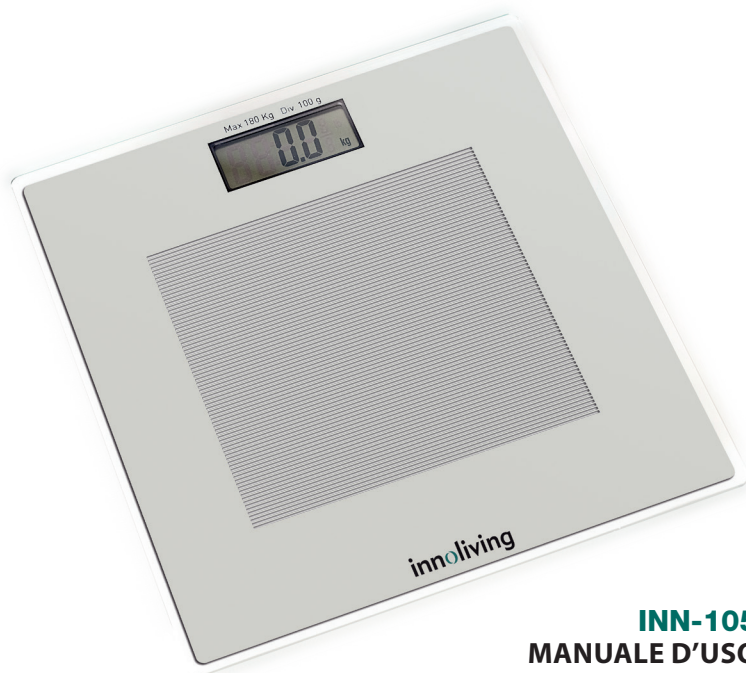INN-105
MANUALE D'USO
USER MANUAL
MANUAL DE USO

MANUALE D'ISTRUZIONE

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INSTRUCTION MANUAL

Thanks a lot for having chosen a INNOLIVING bathroom scale. This is a useful tool in order to check your weight, whenever you want.

TECHNICAL SPECIFICATIONS

Max capacity: **180Kg**

Graduation: 100 g

Tempered glass platform

4 Strain Gauge accuracy sensors

Large and easy to read LCD display

STEP ON: step on the scale to turn it on and the get your weight

Auto switch off

Overloading indicator

Low battery indicator

Kg/lb/st selection

Power supply: 1 lithium battery CR2032 included

OPERATING INSTRUCTION

Open the battery compartment under the scale. Remove the isolating tag and close the battery compartment. Select the weight unit kg/lb/st by the key placed near the battery compartment. Put the scale on a hard, flat surface (avoiding carpets). Step on the scale; the display will show your weight automatically.

The scale will switch off automatically after few seconds.

To repeat the measurement

Step down the scale and wait for the scale turning off; step on the scale again and get your weight.

Attention: Stand still on the scale, distributing your weight evenly. Do not move.

Battery replacement

Open the battery case and insert 1 new battery paying attention to the polarity.

MAINTENANCE AND CARE

Please clean the unit with a damp and soft cloth avoiding any water seepage inside the scale. Do not use chemical detergents. Never deep the product into water. Never step on the scale platform with wet foot. Do not store the scale in an upright position when not in use; this may cause battery exhaustion. Avoid shaking, dropping or shocking the scale. This is a precision instrument and must be handled with care. Never try to repair the scale by yourself but always address to specialised centres. On the contrary, the warranty will be not valid anymore. When the unit is not used for a long period of time, remove the battery, as matter of fact its leaking can damage the appliance. As quality of batteries varies from brand to brand, the product warranty does not include damages to the product caused by battery fluid leakage. When the scale is off, never place or lay any objects on its platform. Always place and store the product away from heating sources, sunlight, dampness and sharp objects.

In case of problems

Check battery have been correctly fitted.

Check that you have selected the weigh unit.

Check whether the scale is on a flat, level floor and not touching against a wall.

Ensure whether you do not step on the scale until the display shows "0.0 Kg".

Low battery indicator

Whether the display shows "Lo" or it does not light up, batteries are flat. Please replace batteries.

Overload indicator

If "EEEE" is displayed, the scale has been overloaded (over 180 kg)

WARNINGS

Use this product only for its intended purpose, as described in this instruction manual such as electronic bathroom scale for domestic use only. Any other use should be considered improper and dangerous. Packaging components (plastic bags, cardboard, polystyrene etc.) must be kept out of the reach of children as they have the potential to be dangerous and their disposal should be carried out according to the current regulations. These scales should not to be used to weigh objects or substances in commercial transactions, to produce medicines, to calculate tolls, tariffs, taxes, premiums, fines, remunerations, indemnities or fees of a similar kind determined by weight. These scales are to be used for keeping one's own personal weight under control, and are not to be used for diagnostic purposes or for medical treatment. Values which may seem unusual or outside the norm should always be discussed with your doctor.

In order to avoid danger to the user of this device, the doses of any drug or therapy prescribed by one's own doctor should in no circumstances be modified on the basis of data calculated by this device.

NB: to avoid possible injury, do not step on the edge of the platform

RATING LABEL

The product is made conforming with all the applicable European regulations.

INFORMATION TO USERS according to Legislative Decree No. 49 of March 14, 2014 "2012/19/ UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE)"

The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014.

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

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