innoliving

**BILANCIA DIGITALE
PESAPERSONA BLUETOOTH**

INN-148
MANUALE D'USO
USER MANUAL
MANUAL DE USO

BILANCIA DIGITALE PESAPERSONA BLUETOOTH

I

MANUALE D'ISTRUZIONE

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INSTRUCTION MANUAL

Thank you for purchasing INNOLIVING digital scale INN-148 with Bluetooth connection. This scale allows you to measure your weight and see the values displayed in your mobile device with iOS or Android system.

Please read the operating instructions carefully and save this manual for future reference.

TECHNICAL SPECIFICATIONS

- Connection with tablets and smartphones with Bluetooth 4.0, Android 4.3 or IOS 8.0 and future versions (Iphone 4S, 5, 5S, 5C, 6, 6 plus, 6S, 6S plus)
- Free download of the App "Innoliving" FROM Apple Store and Google Play.
- Capacity: 180 Kg
- Division: 100 g
- Display of BMI (Body Mass Index)
- Tempered Glass Platform
- 4 StrainGauge Sensors
- Wide and easy to read LCD display
- STEP ON: step on the scale to turn it on and measure your weight.
- Auto switch off
- Overloading Indicator
- Low Battery Indicator
- Kg/lb/st Selection
- Power Supply: 1 lithium battery CR2032 included
- Dimension: 31,9 x 32,8 x 2,2 cm
- Weight: 1.75 Kg

IMPORTANT WARNINGS AND SAFETY RULES

Only use the scales for the purposes indicated in this instruction manual that is as DIGITAL SCALE for domestic use. Any other use is considered improper and therefore dangerous.

Packaging components (plastic bags, cardboard, polystyrene etc.) must be kept out of the reach of children as they are potentially dangerous and their disposal should be carried out according to the current regulations.

These scales should not to be used to weigh objects or substances in commercial transactions, to produce medicines, to calculate tolls, tariffs, taxes, premiums, fines, remunerations, indemnities or fees of a similar kind determined by weight.

These scales are to be used for keeping one's own personal weight under control, and are not to be used for diagnostic purposes or for medical treatment. Values which may seem unusual or outside the norm should always be discussed with your doctor.

In order to avoid danger to the user of this device, the doses of any drug or therapy prescribed by one's own doctor should in no circumstances be modified on the basis of data calculated by this device.

To ensure accurate results, keep the scale away from magnetic fields as this may negatively affect the transfer of data or damage the device.

As a precision instrument, the scale must always be treated with care. Do not drop the scale, do not jump on it or subject it to shocks.

Do not climb on the edge of the scale as this may leads to falls or make measurements inaccurate.

Do not use the scale on tiled or wet floor, as this could lead to falls. Make sure the surface is clean and dry before weighing, as it may become slippery when wet. Do not place the scale on an uneven floor, on a soft surface or on a carpet, because such use could result in unreliable data. Specifications and designs are based on the latest information available at the time of printing and subject to change without notice.

BATTERY INSTALLATION

Open the battery compartment located under the scale. Remove the isolating tag and close the battery compartment.

Rechargeable batteries are not suitable for this type of scale.

When the unit is not used for a long period of time, remove the battery, as matter of fact its leaking can damage the appliance. As quality of batteries varies from brand to brand, the product warranty does not include damages to the product caused by battery acid leakage.

The battery must be inserted respecting the polarity.

UNIT SELECTION

Select the weight unit kg/lb/0:0 lb by the key placed near the battery compartment.

Bluetooth CONNECTION AND INNOLIVING APP

Make sure that the smartphone or tablet has Bluetooth turned on. It is not necessary to go in the settings of Bluetooth to make the synchronization and connect the devices, because the scale is recognized automatically by turning on the App "Innoliving".

Before first use, download and install Innoliving App from the App Store. Launch the App

You must first create a personal profile, by entering all the required data. You can also add your own pictures, clicking on the shape of the profile.

Fig.1

Fig.2

Fig.3

Press **Add Users** on Home to add additional personal profiles, entering all the required data. It is very important that all users of the scale follow the on-screen instructions to register and complete a personal profile as the main data such as height and age are required to measure the BMI (Body Mass Index).

In case you have entered more than one profile, you can choose the user to use by sliding the various profiles in the Home.

To make the measurement, click on the selected profile, or go direct to the Measurement Area.

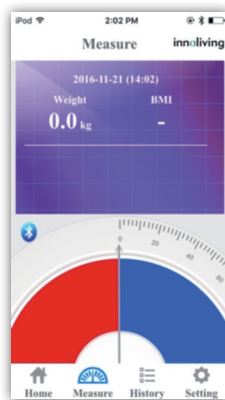
Fig.4

Place the scale on a flat and stable surface (carpets are not recommended). Stand with both feet on the scale for weight measurement, the scale display will show the weight after a few seconds. The scale immediately detects the mobile device and the Bluetooth symbol on the scale lights up and becomes from gray to blue.

Note: to obtain an accurate measurement, when you step on the scale always keep a standing position trying to distribute the weight evenly and avoiding moving. After a few seconds the weight and body mass index BMI will also be displayed on your mobile screen.

Fig.5

The app also displays helpful tips to maintain or improve your physical fitness.

In the "History" area, you can view the measurements of the day, week, month and annual of your own weight. daily performance, weekly, monthly and annual basis of its own weight.

WHAT IS THE BMI?

The Body Mass Index (BMI) is a widely used parameter to get a general assessment of your own body weight.

It correlates the height with the weight of the subject with a simple mathematical formula. It is obtained by dividing the weight in kg of the subject with the square of height in meters. The result of this formula ranks the subject in a specific area of weight: Normal - Underweight - Overweight - Obesity Medium Grade - High Grade obesity.

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ERASING DATA:

You can erase data stored in the user profile, accessing the HISTORY area. Click on the measurement that you want to delete Fig. 6, you will find the details of the detected values fig.7. To delete, press "Delete."

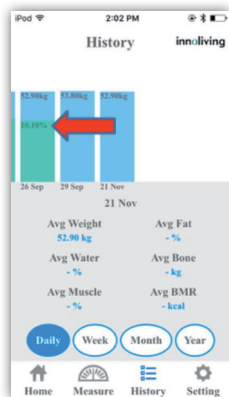
Fig.6

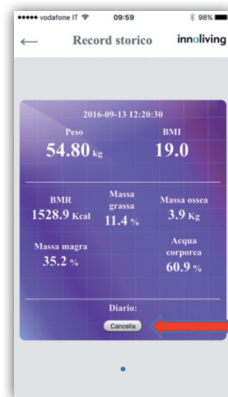
Fig.7

SETTINGS:

In this area you can:

- Change language
- Download and receive an email report
- Synchronize with Apple Health app or Google fit.

Fig.8

To export the data, click "Export Scale Report", the app will display a window where you can specify the period for which you want to receive a measurement report. The report will be sent to the email indicated during registration of the user profile.

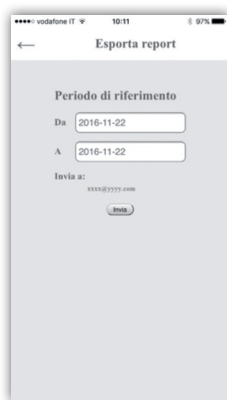
Fig.9

To synchronize the APP with APPLE HEALTH or Google Fit, you must first check that the app HEALTH or Google Fit see among the various sources also Involving, see fig. 3 and 4 here below, and then allow access to the recorded data.

Fig.10

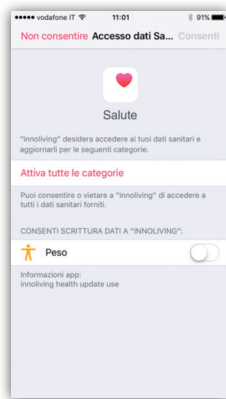
Fig.11

Low battery indicator

If the scale does not turn on or shows “Lo”, the battery may be dead; replace the battery.

Overload indicator

If the display shows “OL”, the scale is overloaded (over 180Kg).

CLEANING AND MAINTENANCE

Clean the scale with a damp cloth to avoid any water infiltration. Never use chemical cleaning agents. Do not immerse in water or other liquids. Do not step on the scale with wet feet. Do not store the scale in an upright position, this fact facilitates the depletion of the batteries. When the balance is not used does not allow anything to rest on the platform weights. Keep away dust and moisture; do not expose to the sun or to excessive temperatures.

Do not open or tamper with the device: repairs must be performed only by authorized personnel; otherwise the warranty expires.

In case of malfunction:

- Check that the batteries are inserted correctly
- Check that the scale is placed on a flat and stable surface and then proceed to the measurement
- Make sure to not climb on the scales before the display will be reset to zero (0.0 Kg)
- If the scale does not connect to the app, make sure the device is connected via Bluetooth to the scale and that the device has the Bluetooth version required by the system.

RATING LABEL:

This product complies with all the applicable European regulations

Read carefully the user manual

INFORMATION TO USERS according to Legislative Decree No. 49 of March 14, 2014
 "2012/19/UE Implementation of the Directive on Waste Electrical and Electronic
 Equipment (WEEE)"

The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014.

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries.

The correct disposal of batteries helps preventing potentially negative consequences on the environment.

The manufacturer Innoliving S.p.a. declares that the BLUETOOTH DIGITAL SCLAE INN-148 complies with Directive 2014/53 / EU. The full text of the EU declaration of conformity is available at the following internet address:

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