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INN
906

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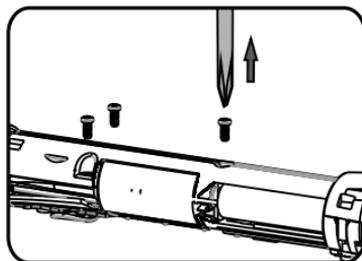
Fig.11

USER MANUAL

PASTA DEL CAPITANO electric toothbrush is a device with oscillating rotary movement, which allows the bristles to completely cover the surface of the teeth, follow the gingival contour and get excellent cleaning results and dental hygiene.

Before starting using the appliance, please read carefully all the instructions below. It is advisable to save this instruction manual for future uses.

IMPORTANT SAFETY INSTRUCTION

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING:

- Before connecting the appliance, check the voltage shown on the rating label corresponds to the local mains voltage.
- This product has been designed to operate with a voltage in the range 100-240V.
- Voltage converters do not guarantee compatibility.
- Do not force the plug to make it fit into the power.
- If the plug does not fit into your outlet, disconnect the product from the power outlet.
- Always disconnect the product from the power outlet after use.
- Never use the product in a state of drowsiness.
- The power cord cannot be replaced. Contact an authorized service center for replacement of the product.
- Do not use this product if the power cord or plug are damaged, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. In case of doubts, please address to a specialised centre.
- If the charging station falls into water, do not try to reach it into the water, but pull out the power plug immediately instead!
- Do not immerse the charger into water or any other liquid.
- Always pull the power cord out of the outlet before cleaning the brush stand or the charging station.
- Keep the unit away from heat, direct sunlight, moisture, sharp objects and similar ones.
- Do not immerse the device in water or other liquid.
- Do not handle with wet hands.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Use only with the charging station and accessories supplied.
- Do not insert any foreign object into any opening on the product.
- Make sure that the power cord does not touch hot surfaces. Do not wrap the power cord around the device.
- Always pulling from the plug and not from the cable.
- Do not dismantle this product at any time other than when disposing of it.
- Do not use outdoors or operate where aerosol products are being used or oxygen is being used or administered. This is a fire hazard.
- Use of this product may interfere with the operation of pacemakers or other internal

medical devices. If you use one of these devices you should consult your physician.

- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Follow the instructions for use, all uses other than described in this manual are considered improper use.

The manufacturer declines all responsibilities for damages to persons or properties resulting from improper use of the device.

This appliance is intended for cleaning teeth and gums. Do not use it for purposes other than those intended.

If you feel pain during use, contact your doctor.

- Keep the device out of reach of children.
- The toothbrush is not suitable for use by children under 3 years of age.
- This device is not a toy. This product is not intended to be used by children without direct adult supervision.


CAUTION: Do not leave the battery compartment open in the presence of children, choking hazard with the ingestion of the batteries.

- The appliance is not intended for use by persons (including children) with mental, physical or sensory capacities reduced, except in cases in which such persons have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Keep all packing material (plastic bags, cardboard, polystyrene, etc..) out of reach of the children, since they are potential sources of danger and must be disposed of as provided by law.
- The product is intended for household use only, do not use the product for commercial purposes.
- Before proceeding with any cleaning or maintenance, disconnect the product from the power outlet.

ATTENTION

- Do not clean the head or the handle in the dishwasher.
- If during the last 2 months you have had interventions on the teeth or gums, ask the advice of your dentist before using the brush.
- Consult your dentist if excessive bleeding after use of the device or if bleeding continues for more than a week.
- In case of oral problems, consult your doctor before using Pasta Del Capitano electric toothbrush.
- Do not use the appliance with crushed or bent bristles. Replace the brush head every two months and whenever appear signs of wear.
- Use only the heads recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the print head with soap and water after each use, to prevent damage to the plastic.

PRODUCT DESCRIPTION:

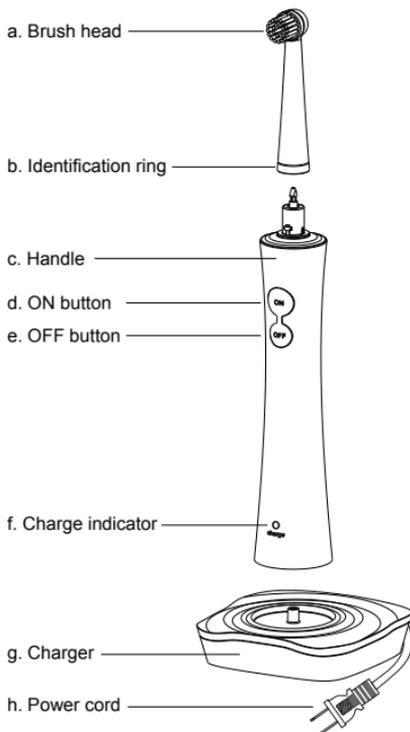HOW TO CHARGE THE TOOTHBRUSH

- Make sure the brush is turned off, before putting it to load.
- Place the handle in the charger (Fig. 1). The charging light will illuminate to indicate the charging process. A full recharge takes about 24 hours.

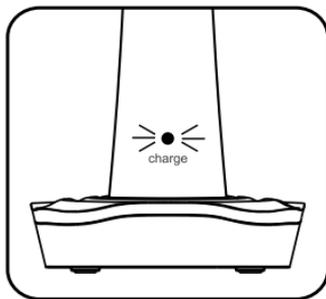
Fig.1

WARNINGS:

- Before using the brush for the first time, you must fully charge the batteries for about 24 hours.
- The toothbrush must be charged 24 hours if:
 - a. you are using the appliance for the first time
 - b. you are not using the unit for more than 3 months
- In order to have the brush fully charged, it is recommended to keep it inserted into the charger when not in use.
- During charging, the toothbrush may become warm, it is normal.
- The useful life of the product decreases with the age of the batteries. If, after fully charging the device, you can use the toothbrush once briefly, it means that the batteries are running low, contact the service center.

HOW TO USE

1. Slide the brush head onto the toothbrush, insert it into the drive shaft catch and turn it clockwise into the drive shaft catch and turn it clockwise a quarter, until it locks into place on handle, to prevent the brush head from becoming detached due to vibration.
2. Moisten the bristles, apply toothpaste.
3. Insert the head attachment into the mouth before switching the unit on.

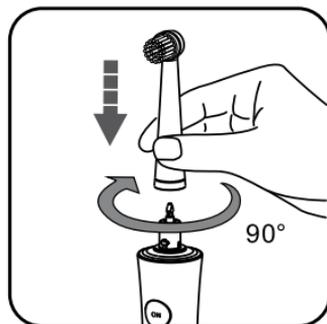
Fig.2

4. Press ON button to turn on the toothbrush (Fig.3).

Brush your teeth following the 4 sections of the dental arches:

outside upper	inside upper
outside lower	inside lower

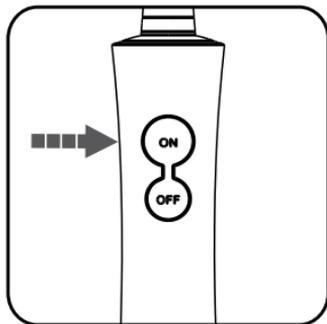
Fig.3

5. Pass the brush slowly on each tooth in each section (Fig. 4), without pressing too hard the brush head.

- Brush each tooth for a few seconds before moving on to the next

- Clean both the teeth, and gums, first the outside then the inside

- Be careful not to brush or press too hard the toothbrush against your teeth and gums

6. Press the OFF button to turn off the brush. Remove the brush head from the mouth only after turning off the device.

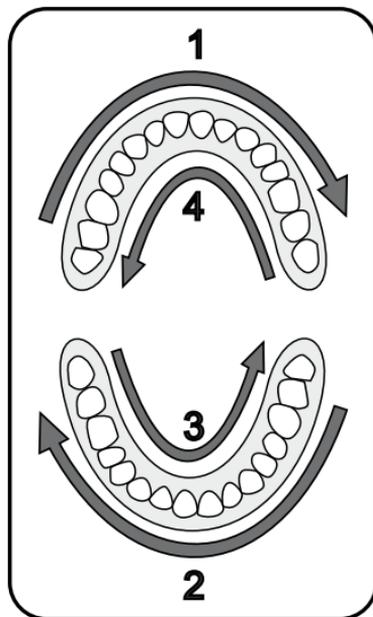
Fig.4

Note:

- In the first days of using an electric toothbrush some users may experience slight bleeding from the gums. This phenomenon should disappear in few days. If bleeding persists, please consult your doctor.
- 2 minutes is considered to be the minimum time needed to brush your teeth properly. For more thorough cleaning, it is advice to prolong the time of brushing a few seconds.

CLEANING AND MAINTENANCE

Thoroughly rinse the brush head after each use. Remove the head from the handle and rinse under lukewarm water, then dry it before you put it back on the handle. The head must be completely dry before you put it in the shaft. Periodically clean the handle and the base as shown below.

1. Remove the head from the handle by turning it anti-clockwise. (fig.5)

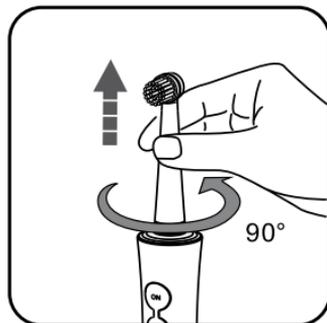
Fig.5

2. Turn on the toothbrush and rinse the vibrant shaft under running water to remove any remaining toothpaste. Turn off the brush.

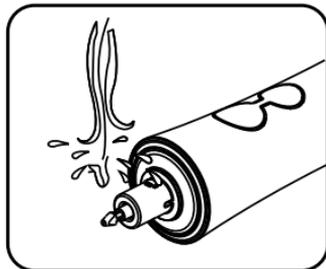
Fig.6

3. Thoroughly rinse the inner part (Fig. 7) and the bristles of the brush head under running water.

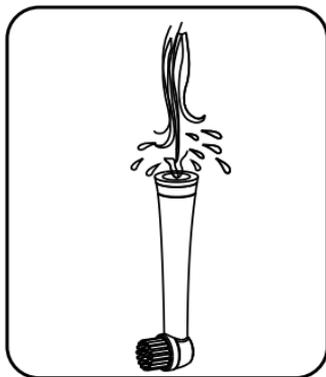
Fig.7

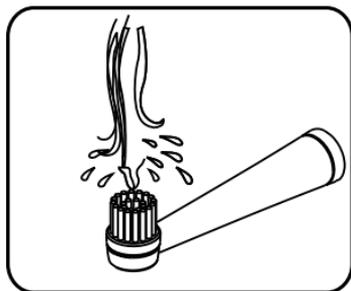
Fig.8

4. Clean the charger base with a damp cloth, then dry. Never immerse the base in water! The base must be completely dry before connecting it to the power outlet. (Fig.9)

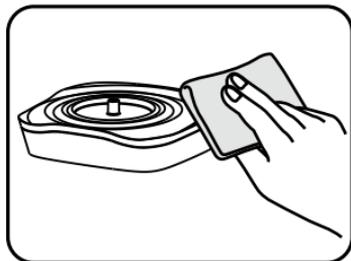
Fig.9

NOTES:

- To keep the brush attachments in an hygienically perfect condition, it is recommended that they are put in a denture cleaning solution approximately every 2-4 weeks in addition to the normal cleaning.
- To obtain a long-lasting, optimal cleaning effect, the brush attachments should be replaced every 2 months or whenever the bristles get damaged, stained or worn folded.
- Any other type of maintenance tasks must be carried out by specialized personnel of an authorized service center.
- If the brush is not used for a long period of time it is recommended to disconnect it from the power outlet and store it in a dry place, away from sources of heat. Keep your toothbrush clean well out of the reach of children. Do not twist the cord around the appliance.

TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush does not work	You have just bought a toothbrush or do not use it for more than 3 months	Load the product for at least 24 hours
Oscillations weak and infrequent	You are pressing the brush head with excessive force on teeth	Gently touch the brush head on the surface of teeth and gums
The toothbrush works only for a few seconds even after charging	The battery has reached the end of its life cycle.	Contact service centre for support
	Charging time was insufficient	Load the product for at least 24 hours

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IPX7

Rev.01_11.2015

Read carefully the instruction manual.

IPX7

Degree of protection from contact with the water.

Device with electrical isolation class II.

This product complies with all the applicable EU directives.

INFORMATION TO USERS according to Legislative Decree No. 49 of March 14, 2014 “ 2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE)”.

The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014.


The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

HOW TO REMOVE THE BATTERY:

The product is powered by 2,4V 600mAh with 2 AAA NiMh batteries included. At the end of the product life, remove the batteries before disposing of the appliance and recycle or dispose of batteries according to local regulations regarding the collection of batteries.

Do not open or disassemble the product, except when you need to dispose of it. Follow the basic safety standards and disconnect the product from the power outlet before removing the batteries.

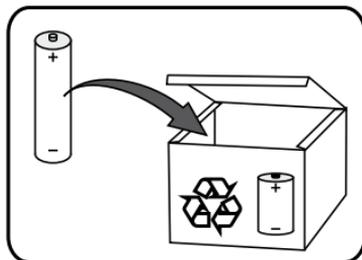Remove the cover from the base of the handle. Use a screwdriver "+", unscrew the 2 screws and open the battery compartment. (Fig.10)

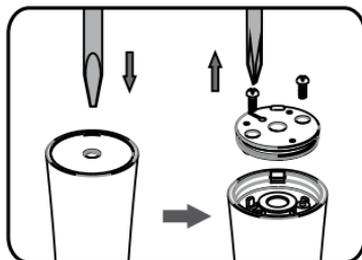
Fig.10

Pry handle handle to remove the battery holder. (Fig.11).

Use a "+" screwdriver to pry the battery out of the battery holder.

CAUTION: Do not connect the appliance to the mains again after you have removed the battery. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

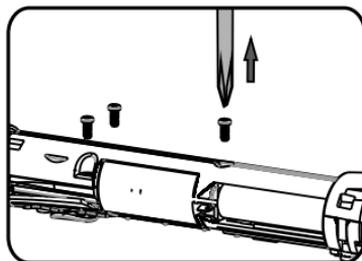
Fig.11

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sotto licenza ufficiale di Farmaceutici Dott. Ciccarelli S.p.A