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RECHARGEABLE ELECTRIC TOOTHBRUSH

GB

### **IMPORTANT SAFETY INSTRUCTIONS**

Read the operating instructions carefully before use and save this user manual for future references.

This device is intended for household use only.

The appliance is not intended to be used by persons (children or adults) with any form of functional disorders, or by persons who do not have sufficient experience or knowledge on how to use it, unless they have received instructions concerning the use of the appliance from someone who is responsible for their safety.

Never immerse the appliance or its base in water or other liquids.

When the charger is connected to the mains do not handle the appliance with wet hands or feet. Only use accessories recommended by the manufacturer.

Do not use the appliance if the main unit or the brush head are damaged.

If more than one person uses the toothbrush, we recommend for hygienic reason that each person use their own brush head.

Do not pull the power cord to remove the plug.

Do not wrap the power cord around the appliance. Do not fold or twist the power cord.

Unplug the appliance in case of a power failure and before proceeding with cleaning operations.

Do not dismantle or attempt to repair the appliance. The repairs should be done by the authorized after sales service or qualified personnel.

This product is not a toy.

Keep the device far away from heating sources.

If the power cord is damaged, it should be replaced by the manufacturer or by the authorized after sales service in order to prevent risks to safety.

This product is not intended for professional use in dental offices. Only use the appliance to clean your teeth and gums; in case of pain or bleeding, contact your physician.

### **TECHNICAL SPECIFICATIONS**

Rated Voltage: 220-240V ~ 50 Hz

Power: 1,6 W

Battery included: Ni-MH 2,4V 600mAh

Charging Time: Approx. 16 hours


Operating Time: 45 minutes

GB

## DESCRIPTION

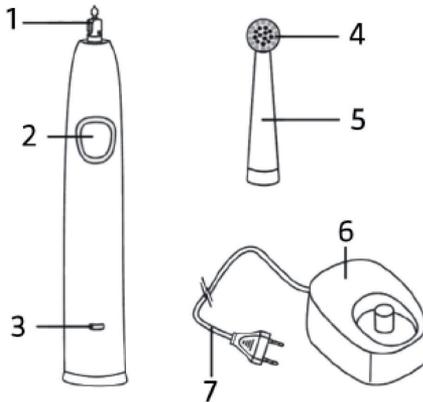1. Holder for Brush Head 2. ON/OFF Button 3. Charge Indicator 4. Bristles 5. Brush Head 6. Charger Unit 7. Power Cord

## HOW TO USE

### Charging the Battery

1. Insert the plug in a power outlet.
2. Place the toothbrush on the charger unit. The charge indicator becomes red.
3. It takes at least 16 hours to complete the charging cycle. The charge indicator becomes green when the battery is fully charged.
4. When the battery is fully charged, you can remove the plug from the power outlet.
5. Charge the battery when you notice a significant decrease in speed, or when the head is completely stopped and the indicator flashes red.


GB

#### NOTE:

- . Turn the toothbrush off during the charging.
- . Charge the battery before using the toothbrush for the first time.
- . Charge the battery when it is discharged, or if the appliance has not been used for more than 90 days.
- . The charger gets warm during the charging. This is quite normal.
- . The battery must be charged with the charger recommended by the manufacturer.

#### Removing the Brush Head

1. Turn the brush head 90 ° to the right (1).
2. Pull off the brush head (2).
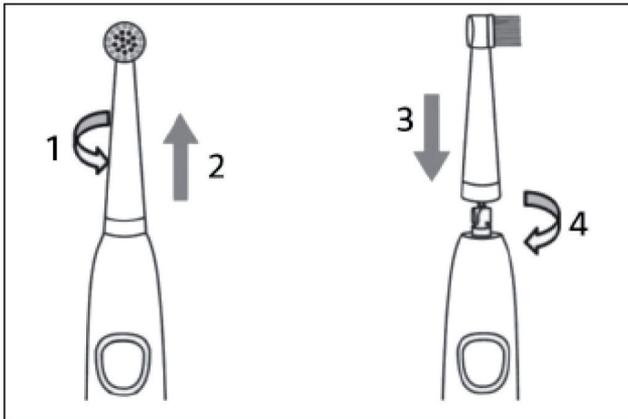
#### Fitting the Brush Head

1. Place the brush head turned 90° to the right in relation to the ON/OFF button.
2. Push the brush head on the shaft (3).
3. Turn the brush head 90° to the left, so that the bristles are facing in the same direction as the ON/OFF button. Check that the brush head is firmly in place (4).

GB

### USE OF THE PRODUCT

1. Put the required amount of toothpaste on the toothbrush.
2. Put the brush head in your mouth and rub the brush against your front teeth to press the toothpaste between the bristles.
3. Press the power switch to turn the toothbrush on and start brushing at the back of your mouth.
4. Brush over all sides of your teeth, between your gums and between your teeth. 2-3 minutes of brushing are usually sufficient.
5. After brushing your teeth, rinse the brush head under running water.

**NOTE:** Do not turn the toothbrush on before putting it in your mouth and switch it off before taking it out to prevent splashing.

Normal toothbrush takes at least 2-3 minutes. The appliance has a timer. After 2 minutes the brush head stops for a moment, and then continues to rotate.

### CLEANING AND MAINTENANCE

After use, rinse the head under running water, remove excess water and store it in an upright position.

Do not immerse the toothbrush or its base in water. Remove the brush head and clean the main unit with a damp cloth.

It is recommended to replace the brush head every 3 months and every time you notice an excessive consumption of the bristles. Use only the spare brush heads model INN-90901 produced with this item.

After use, store the product in a dry place, away from heat sources.

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MADE IN CHINA

GB

Read carefully the instruction manual.

**IPX7**

Degree of protection from contact with water.

Class II insulated electrical appliance

This product is made conforming all the applicable European regulation

**INFORMATION TO USERS** according to Legislative Decree No. 49 of March 14, 2014 “ 2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE) “The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

**innoliving**  
*easy tech, easy life*

Prodotto Da Innoliving Spa  
Via Merloni 9 60131 Ancona  
Tel: 199.188.335 [www.innoliving.it](http://www.innoliving.it)

MADE IN CHINA

Rev.00\_12.2016

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Dott.Ciccarelli S.p.A**