IMPORTANT SAFETY INSTRUCTIONS

- Please read the instructions carefully before use and keep this manual for future reference.
- The appliance is intended for domestic use only.
- This product can be used by children aged from 8 years and above and persons with reduced physical, mental or sensory capabilities, or lack of experience and knowledge, only if they are supervised or have received instructions concerning use of the appliance in a safe way and understand the hazards involved.
- Ensure that children do not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never immerse the appliance or base in water or other liquids.
- When the charger is connected to the mains, do not handle the appliance with wet hands or feet.
- Use only accessories recommended by the manufacturer.
- Do not use the appliance if the body or head is damaged.
- For hygiene reasons, if the product is used by more than one person, it is advisable to use a different head for each user.
- Do not pull the power cord to unplug the product.
- Do not wrap the power cord around the appliance. Do not bend or twist the power cord.
- Unplug the power cord in case of power failure and before cleaning the product.
- Do not attempt to modify or repair the appliance yourself. Repairs must be performed by an authorized service center or qualified personnel.
- This product is not a toy.
- Keep the product away from heat sources.
- The power cord cannot be replaced. If the cord is damaged, the appliance must be scrapped.
- This is not a device for professional use in dental offices. Use it only to clean your teeth and gums. If you experience pain or bleeding, contact your doctor.
- The battery included in the toothbrush is not removable by the user and will be removed by the recovery and treatment center before disposal.
- The toothbrush must be charged using only the charging base provided.


TECHNICAL SPECIFICATIONS

Battery type: Li-ion battery 3.7V/ 600mAh
Nominal voltage: 100-240V ~ 50-60Hz
Power: 3,5W
IP grade: IPX7
Interchangeable heads models: INN-90801NEW, INN-90802, INN-90803, INN-90804, INN-90805, INN-90806
Charging base
LED charge level indicator

FUNCTIONS

- 3 cleaning modes: 1. Clean; 2. Sensitive; 3. Massage.
 - 2 minute timer. The toothbrush automatically turns off after 2 minutes of use.
 - Effective cleaning of lateral and medial teeth, spaces and gingival crevices.
- Suitable for daily use to remove food residues from the oral cavity.

LIST OF COMPONENTS

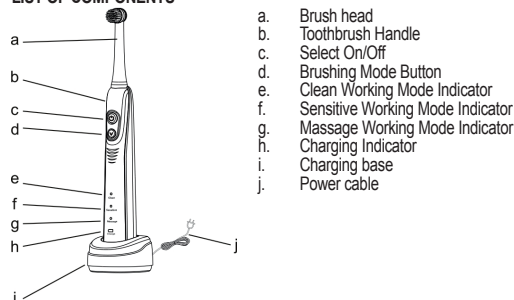CHARGING MODE

MAKE SURE THE SWITCH IS OFF BEFORE CHARGING.

- Place the toothbrush on the charging base, making sure it is properly connected.

- While charging, the blue indicator light flashes continuously.
- Once the charging process is complete, the indicator light remains on.
- If the toothbrush is not used for a long time, the charging indicator may not appear immediately; wait 5-10 minutes for it to return to normal.

Please note that the charging time of the toothbrush is 12 hours.

When the toothbrush head does not rotate strongly, recharge it.

Battery life is 60 minutes in total, equivalent to about 10 days of use with 3 sessions of 2 minutes per day.

Do not charge for more than 24 hours, otherwise it will damage the battery life.

USE

Press the power button.

Select the brushing mode using the mode selection button.

The operating sequence is: clean mode, sensitive mode, and massage mode.

You can select the desired mode even when the device is off using the mode selection button, but the device will turn off if you don't press the power button within 5 seconds.

To switch from clean to sensitive mode, press the brushing mode button.

To use massage mode, press the brushing mode button again.

When you press the brushing mode button, the brushing mode indicator lights up, indicating the selected operating mode.

To turn off the toothbrush, press the power button.

Guide the toothbrush between your teeth to ensure effective cleaning of each tooth.

Hold the brush head in place for a few seconds before moving on to the next tooth.

Brush your teeth starting from the outer surfaces, then the inner surfaces, and finishing with the chewing surfaces. Avoid applying excessive pressure when brushing.

The INN-914 electric toothbrush is equipped with a pressure sensor: if you apply too much pressure (brushing too vigorously), the LED (h) will flash red quickly and you will feel a warning vibration. Reduce the pressure to protect your teeth and gums.

INSTALLING AND REPLACING THE BRUSH HEAD

- Insert the brush head bar directly onto the metal shaft of the handle when installing the new brush head;
- Pull the brush head straight out without rotating it when replacing it with a new brush head;
- Do not install the brush head in the wrong or opposite direction to avoid problems;

CLEANING AND MAINTENANCE

After use, rinse the head under running water, shake off excess water and store the toothbrush upright.

Do not immerse the toothbrush body in water. To clean, remove the head and wipe the toothbrush body with a damp cloth.


We recommend replacing the brush heads every 3 months, or whenever you notice excessive bristle wear. Use only the brush heads designed for this model.

After use, store the product in a dry place, away from heat sources.

RATING LABEL CHARGING BASE

TOOTHBRUSH

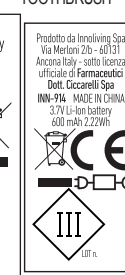Read the instruction manual carefully.

Degree of protection against contact with water.

Class III insulated electrical appliance.

Double insulation

This device complies with all applicable European directives.

INFORMATION FOR USERS pursuant to Legislative Decree No. 49 of 14 March 2014 "Implementation of Directive 2012/19/EU on waste electrical and electronic equipment (WEEE)".

The crossed-out wheeled bin symbol on the appliance indicates that at the end of its useful life the product must be collected separately from other waste.

The user must therefore deliver the equipment, complete with essential components, at the end of its useful life to suitable separate waste collection centers for electronic and electrotechnical waste, or return it to the retailer upon purchase of a new, equivalent type of equipment, on a one-to-one basis, or one-to-zero basis for equipment with a longest side less than 25 cm. Appropriate separate waste collection for subsequent recycling, treatment, and environmentally compatible disposal of the discarded equipment helps avoid possible negative effects on the environment and health and promotes the recycling of the materials from which the equipment is composed. Illegal disposal of the product by the user will result in the application of the administrative penalties set forth in Legislative Decree no. 49 of March 14, 2014.

The battery included in the toothbrush is not removable by the user and will be removed by the recycling center before disposal. The batteries used in this device should be disposed of in the appropriate waste bins at the end of their useful life. Check local regulations regarding battery recycling. Correct disposal of batteries helps avoid negative consequences for the environment and health.

Technical information on the operating modes in accordance with EU Regulation 2023/826:

Mode	Power consumption	Period after which the equipment reaches the described mode
Off	-	-
Standby	0.2309W	When the battery is fully charged and the product is connected to the mains, the device automatically enters standby mode.
Networked Standby	-	-

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